Redi*Stik*TM ECMO Cannulation Trainer


A Collaboration between TCH and Pacific Research Laboratories/Sawbones




Simulation task trainer designed and vetted by frontline innovators at Texas Children's Hospital® in collaboration with Pacific Research Laboratories, Inc./Sawbones®.

RediStik™ ECMO Cannulation Trainer

Anatomically accurate design provides all steps for Pediatric ECMO neck cut down, cannulation, connection to ECMO circuit and de-cannulation in an 18-25 kg patient. Model accommodates 15fr arterial. 17Fr venous and 20 Fr VV cannulas.

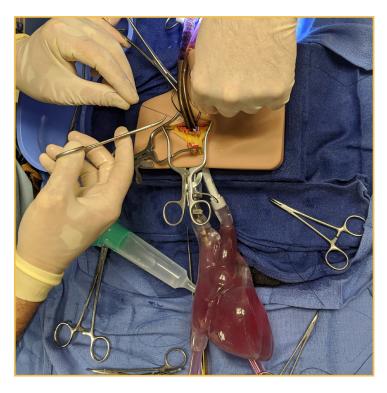

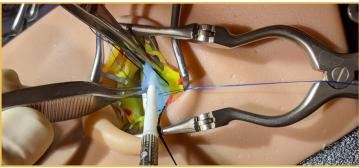



Progress from novice to advanced in ECMO neck cut down, cannulation, connection to ECMO circuit and de-cannulation.







NECK SIMULATOR (SINGLE USE)

- · Right side of neck with anatomically accurate landmarks (mandible, clavicle, SCM).
- · Realistic tissue layers (skin, fat, platysma and SCM) for cutting and spreading to expose the vessels:
 - SCM can be spread, cut, encircled, or laterally retracted to expose vessels.
- · Correct neurovascular (IJV/Carotid/Vagal nerve) positioning, color and feel allows for accurate exposure, cutting, applying traction sutures/tourniquet to control the vessels, and cannula insertion.



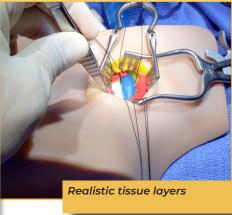

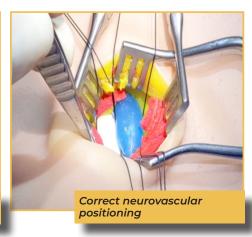


HEART (REUSABLE DESIGN)

- Transparent heart provides visualization during insertion of the wire/cannula.
- Provides closed loop system to connect to the ECMO circuit.
- · An elevated liter bag of saline creates pressure in the system to allow for realistic "bleed back", so a pulsating pump is not necessary.





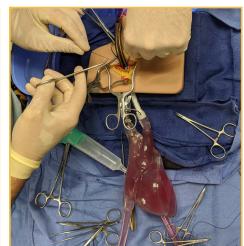



ORDER REFERENCES



RediStik™ ECMO Cannulation Trainer Starter Kit Part #1846-10

Includes:

1 x Duffle bag

1 x Heart/fittings (spare set of fittings for carotid/

IJV + 3 extra stopcocks)

1 x Securement base for heart/neck

1 x Redi*Stik*™ ECMO Heart case

4 x Necks/vessels/vagus nerves (part #1846)


1 x Spare vessel set (one IJV/one carotid — part # 1846-5)

1x Color tabs

1 x Small liquid soap

1x Instructions

(Neck/vessels are for single use. The heart, base for neck/heart, RediStik™ ECMO Heart case, and duffle bag are reusable.)

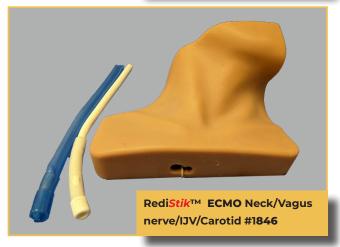

RediStikTM ECMO Neck Replacement (single use) Part #1846

Includes:

1x Neck

1 x Vagus nerve

1 x IJV (Internal Jugular Vein)

1 x Carotid artery

2 x Lubricating jelly pack



RediStikTM ECMO IJV/Carotid Extra Set (single use) Part #1846-5

Includes:

1 x IJV (Internal Jugular Vein)

1 x Carotid artery

2 x Lubricating jelly pack

Visit texaschildrens.org/redistik to access product details, skills video tips and techniques and data outcomes.



in collaboration with

